

New MONTHLY Activity History Report

Reports now listed across the top.

Same Export options available.

EXPORT

PDF/Print

Excel

Raw Export

Choose "Activity By" in New Filters Section (see Page 2).
Also select months and years to view.

Scroll bar will appear if data exceeds page limit.

HH ACTIVITY HISTORY MONTHLY | HH ACTIVITY HISTORY QUARTERLY | HH COMPLIANCE | HH JOB ROLE AND MOMENTS | HH JOB ROLE AND OUTCOME | HH JOB ROLE AND UNIT | **EXPORT**

CUSTOMER: General Healthcare | FACILITY: Demo Hospital OBV In/Out

EDIT FILTERS | SAVE REPORT LINK | REFRESH REPORT

Shifts: All 24 Hours | Locations: All Selected | Activity By: Job Role | Months: All Selected | Job Roles: All Selected | Protocol Moments: All Selected | Years: 1 of 4: 2021

EXPAND LOCATION DETAILS

Job Role	Jan 21		Feb 21		Mar 21		Apr 21		May 21		June 21		July 21		Aug 21		Sept 21		Oct 21		Nov 21		Dec 21
	%	#	%	#	%	#	%	#	%	#	%	#	%	#	%	#	%	#	%	#	%	#	%
Respiratory Therapist	51.8%	575	51.8%	533	54.5%	613	52.0%	548	48.3%	582	49.8%	550	51.7%	580	49.0%	606	48.9%	587	53.8%	513	50.6%	587	48.8%
Physical Therapist	51.8%	583	53.0%	549	50.3%	580	48.5%	582	50.2%	596	53.4%	597	50.4%	570	48.6%	576	48.2%	529	49.9%	565	49.2%	551	49.2%
EVS	50.5%	586	50.0%	528	51.9%	578	47.3%	548	50.8%	587	50.1%	539	49.4%	618	49.3%	588	48.7%	552	53.0%	589	49.3%	552	50.2%
Physician	50.1%	551	50.6%	541	52.0%	586	50.0%	522	49.6%	546	48.6%	553	50.6%	605	51.0%	580	47.6%	534	48.7%	585	53.0%	593	47.5%
Nurse Practitioner	51.9%	603	49.6%	502	44.7%	586	49.4%	539	47.9%	587	48.1%	566	51.1%	560	48.4%	609	49.7%	580	53.3%	546	50.8%	539	48.1%
Aide	52.2%	580	49.8%	512	51.5%	549	50.9%	593	49.2%	563	51.7%	561	53.0%	562	49.1%	550	53.5%	576	44.1%	603	51.9%	572	52.0%
Nurse	51.4%	584	50.8%	526	51.6%	576	49.6%	591	49.6%	585	51.6%	545	52.9%	611	50.0%	586	52.8%	555	47.2%	638	48.7%	557	49.8%
Dietary	51.2%	588	52.8%	508	48.8%	586	48.2%	577	50.6%	603	48.3%	590	48.9%	538	50.2%	562	47.8%	588	49.1%	605	50.6%	549	50.2%
Total	51.4%	4650	51.1%	4199	50.7%	4654	49.5%	4500	49.5%	4649	50.2%	4501	51.0%	4644	49.5%	4657	49.7%	4501	49.7%	4644	50.5%	4500	49.5%

Rows per page: All | 1-8 of 8

New Filtering! Edit Filters now includes all the filtering options in one area. All options are still available for timeframe, job roles and Protocol Moments. Shifts will appear if your facility has shifts set up in the software. All locations can be selected / unselected. as needed.

New "Share Link" feature allows you to save your filtered report by creating a unique URL for future use or to share with other licensed users.

Expand / Hide Location Details allows you to see which floors, units and rooms you have included in the location filters for the report.

Edit Filters

Shifts


Job Roles

Locations

Protocol Moments

Timeframe

CANCEL APPLY

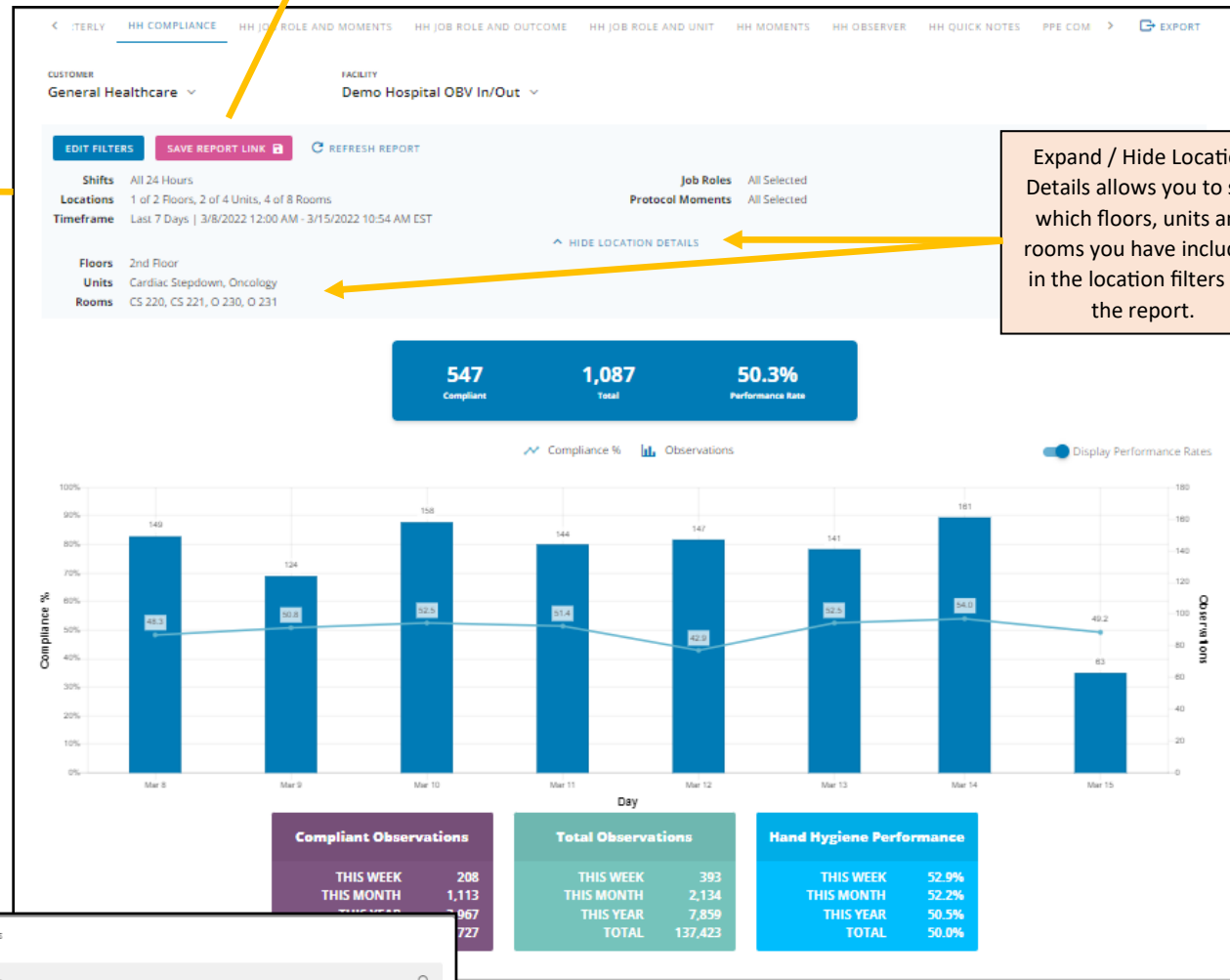

Timeframe

Today

Yesterday

Last 24 Hours

Last 7 Days

Last 30 Days

Last 90 Days

Last 6 Months

Last 12 Months

CANCEL SAVE

Timeframes are the same, and we added a "Yesterday" option.

Locations

Search


Select All Expand All

- 1st Floor SELECT ALL
- CCU SELECT ALL
- ICU SELECT ALL
- 2nd Floor SELECT ALL
- Cardiac Stepdown SELECT ALL
- Oncology SELECT ALL

CANCEL SAVE

To select all the data within a unit, click the "Select All" button next to that unit.

Another noticeable change for some reports is that you no longer have to “opt in” to see the data. For example, when you click on the HH Job Role and Unit Report, all Job Roles are included in the report as a default. You can use the Filters to deselect any job roles you’d like to exclude.

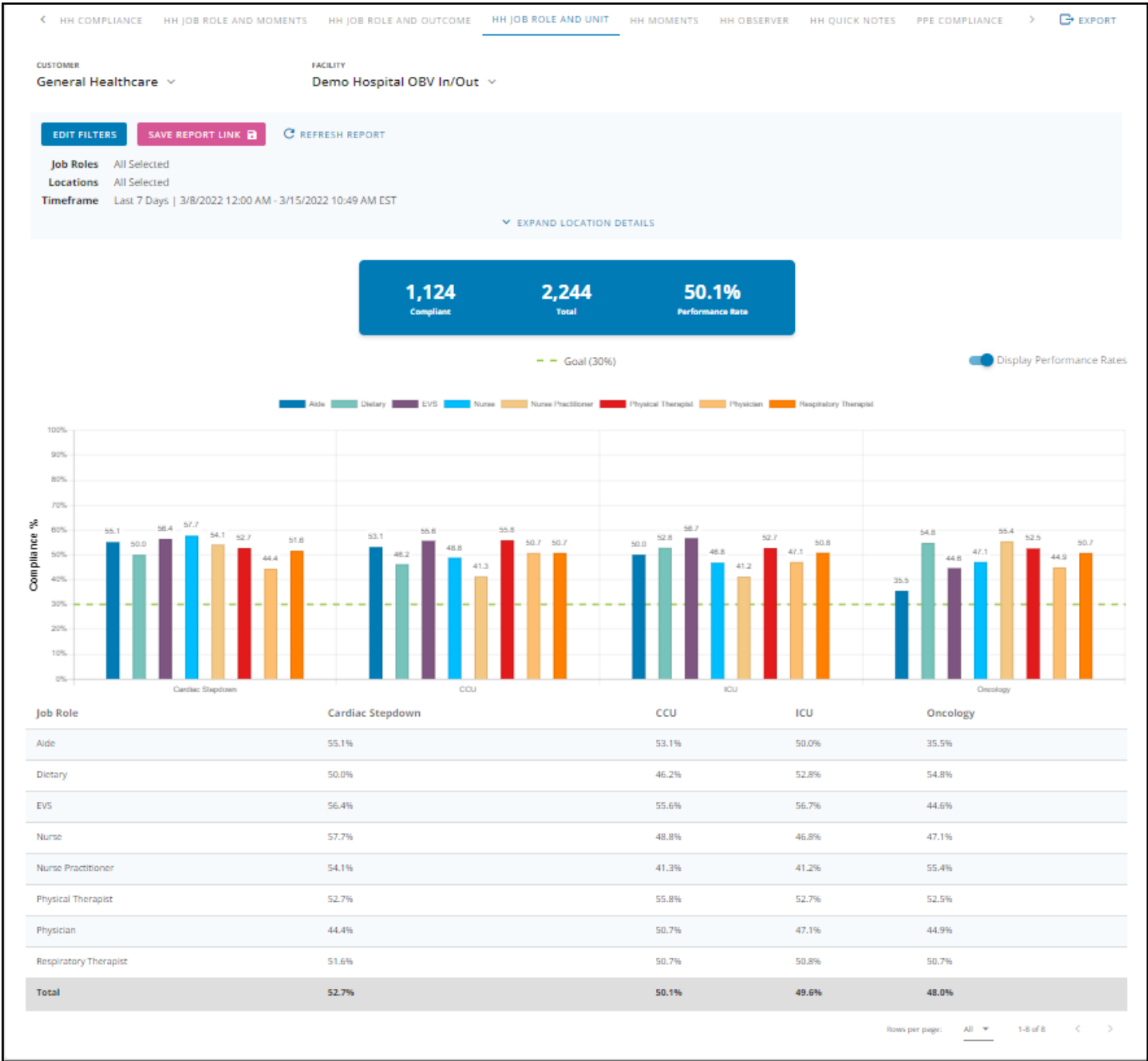

Need Help? Email us at smartlinksupport@gojo.com!